


# FINDER

Finder measures the COSMETIC LINE, a linear measurement from the biting edge of the upper front teeth to the tip of the nose. The norm for the cosmetic line is the age of the child plus 23. This millimeter measurement varies only slightly with body size. An adult female should have a 38 to 42 mm. measurement, and an adult male should measure from 40 to 44 mm.

The following is a table of ideal COSMETIC LINE measurements for a given age. An increase of 5 or more millimeters from these standards at any given age is an indicator of undesirable facial growth.

AGE	MALE	FEMALE
5	28mm.	28mm.
6	29mm.	29mm.
7	30mm.	30mm.
8	31mm.	31mm.
9	32mm.	32mm.
10	33mm.	33mm.
11	34mm.	34mm.
12	35mm.	35mm.
13	36mm.	36mm.
14	37mm.	37mm.
15	38mm.	38mm.
16	39mm.	Adult female 38-42 min
17	40mm.	
Adult male 40-44 min		

1. Hook the measuring device on the upper front teeth.
2. Place the arm to the tip of the nose.
3. Read the mm. measurement at the most forward point of the nose.
4. Compare the mm. measurement and age of the child to the table above.

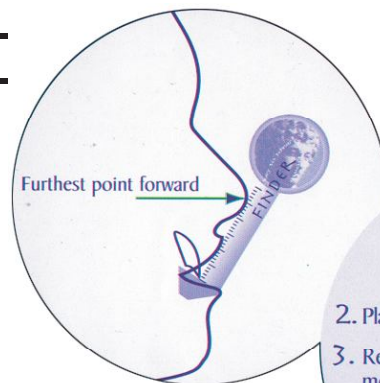

1. Hook the measuring device on the upper front teeth.
2. Place the arm to the tip of the nose.
3. Read the mm. measurement at the most forward point of the nose.
4. Compare the mm. measurement and age of the child to the table above.

